

DIGITAL ECG SOLUTION BY MIDMARK



IQ
ecg





The IQecg™ from Midmark is an easy-to-use, digital 12-lead resting ECG solution that helps you quickly administer testing and access patient information you need to make informed cardiac assessments and provide more efficient patient care.

INNOVATIVE FEATURES HELP IMPROVE PATIENT CARE

Efficiency

- Lightweight, portable design, weighing less than 11 oz.
- View and monitor live, 12-lead ECG in selectable formats on screen
- Capture abnormal ECG as it appears at the touch of a button
- Midmark 12-lead Resting Electrocardiogram Analysis Program automatically analyzes and interprets ECG data for neonate, pediatric and adult patients
- Review and edit ECG results on-screen with free text or customizable quick edit statements
- Detects possible arm lead reversal during the test so the technician can take appropriate action while the patient is still connected

Technology

- Save time and reduce data entry errors by populating patient information directly from your EHR/EMR software
- Facilitate a paperless office strategy and eliminate the costs of thermal paper and the hassles of scanning and filing patient charts
- Stores binary ECG waveforms and displays in digital data instead of static images, resulting in smaller file sizes and providing accurate viewing and measurements
- Quickly assess ECG changes with direct overlay serial comparison of up to four ECG tests
- Enhance your diagnosis with electronic tools like zoom and on-screen calipers

Adaptability

- Functions as a standalone product with our IQmanager™ software or seamlessly integrates with more than 80 of today's top EHR/EMR solutions
- Integration allows you to perform tests within the EHR/EMR and store the information in the patient's electronic chart
- IQpath™ software solution offers uninterrupted testing, reduced latency issues and smooth data flow for facilities operating in a thin-client environment
- Part of a family of products created to grow with your practice



Lightweight, portable design, weighing less than 11 oz.

IQecg™: An ECG that grows with your practice

Functions as a standalone product with our IQmanager™ software or seamlessly integrates with more than 80 of today's top EHR/EMR solutions.




An easy-to-use, 12-lead digital ECG that helps to quickly administer testing.



Review and edit ECG results on-screen.



Lightweight, portable design makes moving from room to room quick and easy.

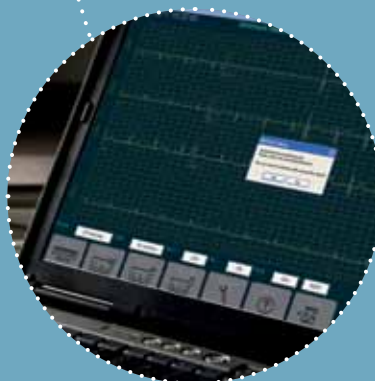

Detects and warns possible arm lead reversal while in live ECG screen.

Make your practice more efficient with an IQ Digital Diagnostic System

Workflow efficiency and improvements in patient care are just a few of the benefits that can result from implementing an IQ Digital Diagnostic System. Call us today and we'll put together a customized suite of solutions to grow as your practice grows.

Diagnostic Products				
	EHR/EMR Connectivity	Single Platform (IQmanager™)	Standalone	Thin - Client (IQpath™)
IQecg™	■	■		■
IQholter™	■	■		■
IQholter™EX	■	■		■
IQholter™EP	■	■		■
IQstress™	■	■		
IQspiro™	■	■		■
IQvitals™	■	■	■	■
IQvitals PC™	■	■		■

Minimum System Requirements

- Windows® Vista, Windows® XP , Windows® 7, (32-bit and 64-bit)
- Intel® Core™ 2 Duo Processor E4300 or 64-bit processor or faster
- 2 GB of RAM
- 2 GB available hard disk space
- USB port
- Mouse, keyboard and high resolution display monitor

Contact us or visit midmark.com for more detailed system requirements.

Windows is a registered trademark of Microsoft Corporation.

Intel Core is a trademark of Intel Corporation in the U.S. and/or other countries.

Models


4-000-0061 IQecg™ Digital ECG USB Version

4-000-0060 IQecg™ Digital ECG Serial Version

Accessories



Universal ECG Clips (3mm, 4mm, & Snap)

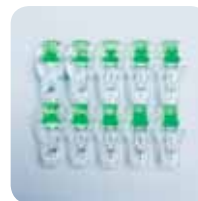

Clear ECG Clips (4mm & Snap)



Disposable EZ-Trode ECG Electrodes



IQmanager™ Software



IQcart™



RETHINK THE CLINICAL SPACE

At the core of Midmark is a new way to think about caregiver/patient interaction and the delivery of efficient patient care. We are dedicated to working together with you to redefine the future of the clinical space.

midmarkclinicalsolutions.com

For information, pricing or to place an order, contact us at:

Medical Device Depot
Inc.

www.medicaldevicedepot.com
1-877-646-3300