

The PowerLifter I and II provide solutions to your organization's elimination of workplace hazards and ergonomic controls. These unique products are specifically designed to neutralize irrigation bag weight and lessen the burden of lifting heavy solutions during surgeries. Used in hospitals, clinics and surgery centers across the United States and internationally during arthroscopic, cystoscopic, orthopedic and TUR surgeries and procedures. These units are simple to use and can be easily operated by one person.

#BEAM

Power Lifter®


- Safe air-operated system featuring an internal, rechargeable 12v battery
- Easy to operate, simply depress the "up" foot pedal to elevate as much as 40 lbs to a height of 9' in seconds
- Depress the "down" foot pedal and the pole lowers quickly and smoothly, making solution changeovers a snap
- Heavy duty base provides constant stability, at maximum and minimum capacity
- Holds up to four 5000 cc irrigation containers
- Battery operated convenience eliminates hazardous electric cords on O.R. floor



The stainless hanger assembly accommodates four irrigation containers up to 5000 cc each. Easy on and off hooks makes solution change easy and quick.



The aluminum locking collar is designed for strength to secure the pole in its raised position. An inner sleeve protects the inner tube against marring.

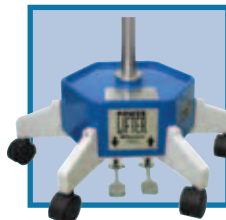

The six nylon nonconductive dual wheel casters support the base and allow smooth transport with minimal effort. The base features a reservoir to collect fluid and reduce spillage on the O.R. floor.



The external battery charger plugs into the access port on the base of the POWER LIFTER®. (Note: The POWER LIFTER® will not work while the unit is plugged in. We suggest that a recharge schedule be established to prolong the battery life.)



Medical Device Depot inc.

www.medicaldevicedepot.com
1-877-646-3300

POWER LIFTER.
#741314

Specifications:

Capacity - Up to four 5000 cc bags/containers
Material - Anodized aluminum poles. Stainless steel hanger assembly. Cast aluminum base. 6 non conductive casters
Lift Mechanism - Pneumatic
Height - Min. 5' 7" from floor to top of hanger - Max. 9' from floor to top of hanger
Shipping Weight - Base - 48 lbs Pole assembly - 8 lbs
Operating Voltage - 12V DC (nom.)
Recharge Power - 115V 60HZ, approx. 180MA

	Power Lifter® Irrigation Stand #741314	Power Lifter® II Irrigation Stand #741315	Low Voltage Light #741316
Description	57 lbs.	78 lbs.	Factory Attached
Price/Each			