

Precise information.
Easier diagnoses. Better patient care.



INTRODUCING

The Welch Allyn SpiroPerfect™ PC-Based Spirometer

Medical Device Depot
Inc.

www.medicaldevicedepot.com
1-877-646-3300

WelchAllyn®

Advancing Frontline Care™

Accurate, integrated patient information at your fingertips

In seconds, you and your staff can see and compare accurate results from pulmonary function tests with the friendly, featherweight Welch Allyn SpiroPerfect™ PC-Based Spirometer. You can clearly evaluate pulmonary obstruction or restriction, increase success rates of patients who are undergoing smoking cessation, show progress, and help pediatric patients comply with test procedures.

- Fast, reliable pre- and post-bronchodilator FVC, SVC and MVV testing with automatic interpretation and reversibility analysis
 - > View real-time Flow/Volume Loop and Volume/Time graphs
 - > Instantly verify quality and variability of proper test performance
 - > Motivate children with incentive graphic—for pediatric patient coaching
 - > Reduce risk of cross contamination with single-use disposable transducers
 - > Select single and multiple flow calibration protocols
- Easily monitor, review, store and print on your existing computer with the Welch Allyn CardioPerfect™ Workstation
 - > Choose from multiple predicted norms and report formats
 - > Compare previous tests with the Trending tool
 - > Make available to multiple staff via computer network
 - > Interface with Electronic Medical Records or patient management systems for onsite or offsite review of key data. Telemedicine option is available for e-mail transfer.


Medical Device Depot
Inc.
www.medicaldevicedepot.com
1-877-646-3300

Integrated. Customizable. Intuitive.

The screenshot displays the Welch Allyn CardioPerfect software interface. On the left, a 'Patient records' list shows various patient names and IDs. The main area features two graphs: 'Flow/Volume Loop' and 'Selected Flow/Volume'. Below the graphs is an 'Interpretation' section with text and a 'Comment' field. To the right, a 'Customizable parameters' table is visible, showing various lung function parameters and their values across different test attempts.

Effort	Units	NHANES III	Post 3	Post 2	Post 1	Pre 3	Pre 2	Pre 1	Final
FVC	L	3.0	2.0	2.0	2.8	3.6	3.4	3.2	3.6
FEV1	L	2.3	2.8	2.7	2.5	3.2	3.1	2.9	3.2
FEV1%	%	78%	92%	94%	92%	90%	91%	96%	90%
FEV6	L	2.1	3.0	2.8	2.8	3.6	3.4	3.2	3.6
PEF	L/s	5.8	6.8	7.0	7.2	7.9	7.0	7.1	7.9
FEF25-75	L/s		3.6	3.6	3.4	4.2	4.1	3.7	4.2

- 1 Store, retrieve, print or edit **patient records** using the Welch Allyn CardioPerfect Workstation database.
- 2 Easily **identify** all diagnostic test types and dates each procedure was performed.
- 3 **Automatic interpretation** follows ATS Interpretative Strategies and test reproducibility criteria.
- 4 See a quick snapshot of patient results. Select from 6 parameters to show real-time testing and summary. View measurement tables with 19 fully **customizable parameters**.
- 5 **Display real time curves** as Flow/Volume, Volume/Time, Flow/Time or Pre, Post. All efforts or best effort can be displayed on the same graph for easy comparison.

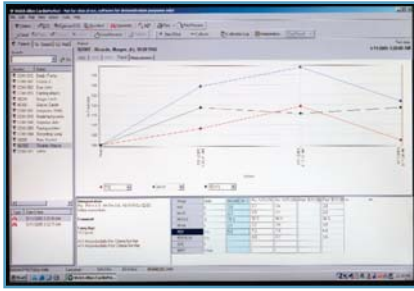

Compare tests done on different days and clearly see progress. View three parameters at once as a percentage of predictive values or of absolute value.



Enjoy the convenience of our accurate patented **Disposable Flow Transducers**. Designed for use by an adult or child, these single-use disposables minimize the risk of cross contamination. No sterilization, no cleanup—just safety and savings.



Coach, motivate and make testing more fun. Because children understand how to help the fireman put out the fire, the **Incentive Screen** helps children understand your goals and comply with testing protocols.

Get connected to your EMR through Welch Allyn's CardioPerfect™ Workstation Software



Welch Allyn SpiroPerfect™
PC-Based Spirometer



Welch Allyn CardioPerfect™
Workstation Software



Transfers Test Data
to Electronic Medical
Records Systems

Medical Device Depot
Inc.
www.medicaldevicedepot.com
1-877-646-3300

The Welch Allyn SpiroPerfect™ PC-Based Spirometer downloads information to the Welch Allyn CardioPerfect™ Workstation Software, which provides paperless efficiency as well as connectivity to many leading EMR systems. This software puts patient records from a variety of Welch Allyn diagnostic instruments at your fingertips, giving you the ability to store, access, review, compare and edit diagnostic test data from any computer on your network.

- Retrieve records with the touch of a key, and review and edit right on the screen
- Store thousands of records with ease
- Access information from your exam rooms, your office, or at home
- Fax or e-mail reports directly for consults
- Add additional cardiopulmonary diagnostic testing devices like ECG, ABPM, or Holter as your practice grows
- Copy and paste spiromographs into Microsoft® Word
- Integrate with other patient information systems
- Uses passwords to assign/limit tasks so that data integrity is maintained

Welch Allyn SpiroPerfect™ PC-Based Spirometer

Technical Specifications

SpiroPerfect	Computer-based full diagnostic spirometer	Quality Checks	Effort Acceptability and Test Reproducibility checks
Tests	FVC, SVC, MVV, Pre- and Post-BD		Audio and visual incentive for assistance in coaching patients
Sensor Type	Pneumotach	Connectivity	Part of the CardioPerfect™ Diagnostic Suite
Power Equipment	None, obtained from USB port		Export compatible with most electronic medical records programs
Accuracy	Meets or exceeds 1994 ATS standard		Available in multiuser network
Reproducibility	Meets or exceeds 1994 ATS standard		Telemedicine option for e-mail transfer
Volume Range	0-14 L		
Flow Range	±14 L/sec		
Software	Microsoft® Windows® 2000 and XP compatible	Minimum Computer Requirements (Recommended)	
Predictive Norms	Berglund 1963, Crapo 1981, Composite, Dockery 1983, ECCS/Quanjer 1993, Falaschetti 1995, Gulsvik 2001, Hedenström 1986, Hsu 1979, Knudson 1976, Koillinen 1998, Kory 1961, Morris 1971, NHANES III 1999 [Default], Polgar 1971, Schoenberg 1978, Solymar 1980, Viljanen 1981, Zapletal 1969 <small>Additional predictive norms can be added upon customer request</small>	Processor	Intel® Pentium® III (Pentium IV or higher)
Interpretation	1991 ATS Interpretation Standards. Automatic can be disabled. Manual interpretation available. Lung age calculation	Hard Drive	2 GB (4 GB or more)
Reports	Multiple report templates including FVC, SVC, MVV, measurement, and trend	RAM	256 MB (512 MB or higher)
Incentive Graphic Parameters	Fireman FVC, FIVC, FIV1, FIV1%, FEV0.5, FEV1, FEV2, FEV3, FEV5, FEV6, FEV0.5%, FEV1%, FEV2%, FEV3%, FEV5%, FEV6%, PEF, FEF25 (MEF75), FEF50 (MEF50), FEF75 (MEF25), FEF0.2-1.2, FEF25-75, FEF75-85, PIF, FIF50, FEF50/FIF50, FET SVC, ERV, IRV, VT, IC, BF, MV, Tin, Tex, Tin/Tex	Monitor Resolution	(1024 x 768 true color)
		Operating Systems	Microsoft® Windows® 2000 or XP
		CD-ROM	Required for installation
		Printers	Color printer recommended. Compatible with inkjet and laser printers. Laser printers require a minimum of 2 MB of memory
		Performance Standards	
			SpiroPerfect meets or exceeds the following standards:
			ATS—American Thoracic Society, Standardization of Spirometry, 1994 Update
			ERS—European Respiratory Society
			NIOSH—National Institute for Occupational Safety and Health
			OSHA—Occupational Safety and Health Administration
			SSD—Social Security Disability
			ECCS—European Community for Coal & Steel

Ordering Information

SPIRO-S	SpiroPerfect™ PC-Based Spirometer Module: PC Spirometer Module with 3 L calibration syringe. Includes: sensor, pressure tubing, 4 flow transducers, electronic operator's manual, Welch Allyn CardioPerfect™ Workstation Software CD, 2 training posters, and quick reference card
SPIRO-NS	SpiroPerfect™ PC-Based Spirometer Module: PC Spirometer Module without 3 L calibration syringe. Includes sensor, pressure tubing, 4 flow transducers, electronic operator's manual, Welch Allyn CardioPerfect™ Workstation Software CD, 2 training posters, and quick reference card
703418	Disposable Spirometer Flow Transducer, 25 pk
703419	Disposable Spirometer Flow Transducer, 100 pk
703415	Pressure Tubing, 2 m
100680	Nose Clip (25 pk)
703480	Calibration Syringe (3 liters)

CPT Codes

94010	Basic Spirometry—Spirometry including graphic record, total and timed vital capacity, expiratory flow rate measurement(s), with or without maximal voluntary ventilation
94060	Pre- and Post-BD—Bronchospasm evaluation: spirometry as in 94010, before and after bronchodilator (aerosol or parental), or before and after exercise

To learn more about the Welch Allyn SpiroPerfect PC-Based Spirometer or to request a demonstration, contact your Welch Allyn representative. You'll find a solution that not only improves patient care, but provides priceless peace of mind.

Medical Device Depot
Inc.
www.medicaldevicedepot.com
1-877-646-3300