




CARDIOVIT MS-2010

Defining the new dimensions of electrocardiography



Discover our revolutionary new electrocardiograph with Full Touch Screen Technology.

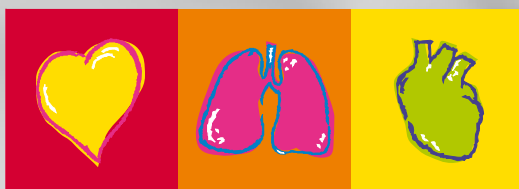

SCHILLER

The Art of Diagnostics

CARDIOVIT MS-2010

See and touch the quality

See the SCHILLER MS-2010 and get a clear sense of the future of cardiography – uniting precision performance and sophisticated ergonomics.

Prepare to be wowed at your first encounter and for a similar sensation every time you use it. MS-2010 is designed for large volumes of work – in hospitals where the ECG volume is high and speed and quality are paramount.

Simply touch the large (10.4 in) high-resolution color display and you're in business – recording, selecting and printing the highest quality ECGs faster than you can imagine. Since MS-2010 exports ECGs in XML, PDF and DICOM format, you can transmit recordings to any system that supports this industry standard. From accurate and complete data entries to comprehensive data fields that prompt users to collect required information, every detail of MS-2010 is designed to streamline ECG testing. And when you're responsible for hundreds of ECGs across miles of hospital terrain every day, attention to detail makes all the difference.

SCHILLER AG
pioneering company
since 1974

Through an ingenious un
a whole new loc



Advanced processing... no waiting

- Instant start-up provides immediate availability
- ECG processed, interpreted and printed in less than 3 seconds after acquisition

Synchronize time

Auto set the MS-2010 time with your hospital time master to obtain accurate documentation of your patient's clinical history.

Fast entry

Enter patient data quickly and accurately as you choose: type on touch screen keyboard or scan the bar code. You can also download patient data and orders from a worklist via SEMA or HIS.

Clinical excellence

Get complete clinical information with the SCHILLER Interpretation software with **thrombolysis** software or with advanced ECG tools such as **SAECG** or **HRV**.

Never miss a beat

Capture and store 5 minutes of uninterrupted 12-lead data to document intermittent arrhythmias.

High-volume storage

- Storage of up to 350 ECGs
- Full disclosure storage of up to 5 minutes for easy acquisition of clean and stable ECG traces .

ion of form and function, the SCHILLER MS-2010 delivers
ok, feel and standard for electrocardiographs.



Display

10.4" full color touch

Network

One of the most advanced networking support found in any electrocardiograph.

Speed

Instant-on and ECG printout in less than 3 seconds

DICOM

First high-end electrocardiograph to support DICOM standard



One-touch workflow

MS-2010 can be configured to automatically print, save, transfer, and retrieve a previous ECG - at the touch of a button.

Connect:

- via LAN, WLAN or GPRS to SEMA or HIS.
- PDF export of resting ECG reports

Automated workflow:

- Acquire, print, save, transfer and retrieve diagnostic ECGs to and from SEMA or HIS - all automatically
- Support for HL7 worklists

Stay connected

Built on a s standard platform, the MS-2010 fits seamlessly into your existing IT infrastructure - so that you are always connected.

Guaranteed safety

The MS-2010 offers you safe wireless connections via standard WLAN protocols such as 802.11 g, WPA, WPA II or WEP. WLAN communication is also available with invisible SSID.

SCHILLER devices support DHCP or static IP to protect your patient's privacy. Even if all safety measures should fail during data transmission, the patient data cannot be read because it is encrypted by SCHILLER.

