

seca 874

Flat scale for mobile use



- Mother/child function for trouble-free weighing of small children.
- Large leveling bases for increased stability.
- Low device weight and battery power for mobile use.
- Double display for readout from different positions.
- Two push buttons for effortless operation with the foot.

seca 874:

The scale that thinks of people.

The seca 874 can do everything that a scale in mobile use should do and then some. It features extremely robust construction and easy-to-use design with long-life batteries and precision unaffected by humidity levels. The scale is so simple that you need only the tip of your toes to operate it.



Great stability.

The seca 874 is very light at only 8.9 lbs, but also very tough. Reinforced by steel elements in the frame, the scale also boasts a top-quality ribbed surface, extra-strength two-component ABS plastic and Plexiglas. Small children and adults up to 440 lbs can be weighed accurately. Especially large leveling bases keep the scale stable on any type of flooring.


Simple operation.

When conducting one examination after the other, there's no need to bend over repeatedly, thanks to the two push buttons on the seca 874 that can be turned on and off with a toe tap. Yet another labor-saving feature is the mother/child function that weighs a child held in its mother's arms. Transport of the scale is made easy by the integrated handle, very low device weight of 7.5 lbs and a carrier bag of feather-light nylon, which is included in the supply schedule.

Intelligent details.

Special mention should be made of the double display that lets the patient and medical personnel read the results from two different directions at the same time. Thanks to a protected toggle switch and a manual power switch, which prevents unintended activation during transport, the function of the scale is always ensured even under difficult conditions.

Both doctor and patient can see the result at the same time in the double display.

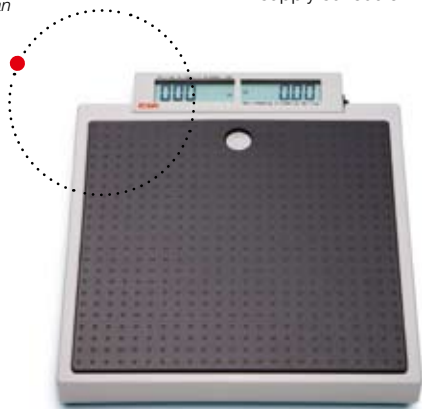

The fit at design and slip-proof surface ensure easy access and a secure stance.



The large leveling bases on the bottom ensure stability on any type of flooring.



The two push buttons on the rear side of the scale can be operated with the foot.

seca 874

Technical data

- Capacity: 440 lbs / 200 kg
- Graduation: 0.1 lbs < 330 lbs > 0.2 lbs / 50 g < 150 kg > 100 g
- Dimensions (WxHxD): 12.6 x 2.4 x 14.2" / 321 x 60 x 361 mm
- Device weight: 8.9 lbs / 4 kg
- Power supply: batteries
- Functions: mother/child-function, auto-HOLD, lbs/kg switch-over, automatic switch-off
- Optional: stadiometer seca 217 combined with adapter element seca 437, carrying case seca 414 (for combination with stadiometer and adapter element) and seca 421

Medical Device Depot Inc.

www.medicaldevicedepot.com

1-877-646-3300