

Power Lifter®

The Power Lifter® I and II provide healthcare workers with hazard-free, ergonomic lifting of solutions during surgical procedures. These irrigation stands neutralize bag weight to ease raising and lowering and prevent workplace hazards. Great for Arthroscopic, Cystoscopic, Orthopedic and TUR surgeries and procedures! These units are simple to use and can be easily operated by one person.


Power Lifter® - Single Column System

- Safe air-operated system featuring an internal, rechargeable 12V battery

Easy to operate:

- Depress the “UP” foot pedal to raise 5,000 cc per hook (total 40 lbs) to a height of 9’
- Depress the “DOWN” foot pedal to smoothly lower solution containers and make changeovers
- Heavy cast aluminum base provides stability at maximum and minimum capacity
- 4-hook independent hanger assembly safely holds a total of 20,000 cc irrigation containers
- Battery operated convenience eliminates hazardous electric cords on O.R. floor

Product	Description	Weight	Height
#741314	... Power Lifter® Irrigation Stand	57 lbs.	Max 9’
#741316	... Low Voltage Light Kit		

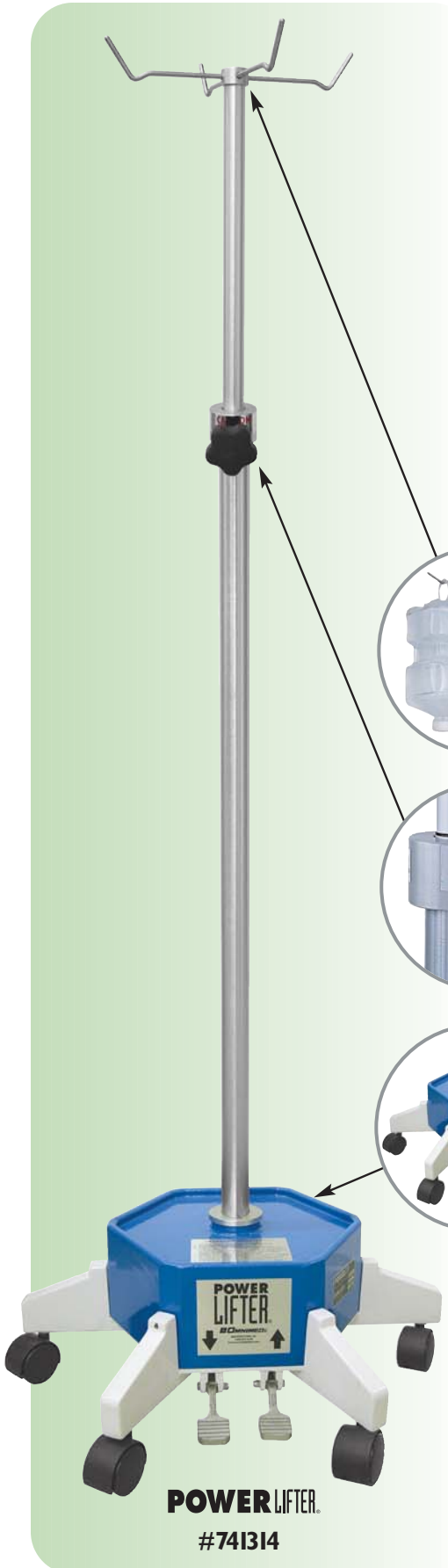

POWER LIFTER.
#741314

The stainless hanger assembly accommodates four irrigation containers up to 5000 cc each. Easy on and off hooks make solution change easy and quick.

The aluminum locking collar is designed for strength to secure the pole in its raised position. An inner sleeve protects the inner tube against marring.

The six nylon non-conductive dual wheel casters support the base and allow smooth transport with minimal effort. The base features a reservoir to collect fluid and reduce spillage on the O.R. floor.

Specifications for Power Lifter® I & II:

- Capacity** - Up to four 5000 cc containers (one bag per hook)
- Material** - Anodized aluminum poles / Stainless steel hanger assembly / Cast aluminum base 6 non-conductive casters
- Lift Mechanism** - Pneumatic
- Height** - Min. 5’ 7” to Max. 9’ from floor to top of hanger
- Shipping Weight** - Power Lifter® / Base - 48 lbs Pole assembly - 8 lbs
Power Lifter® II / Base - 48 lbs Pole assembly - 32 lbs
- Operating Voltage** - 12V DC (nom.)
- Recharge Power** - 115V 60Hz, approx. 180MA

Power Lifter II®

Power Lifter® II - Dual Column System

- Safe air-operated system featuring an internal, rechargeable 12V battery

Easy to operate:

- Depress the “UP” foot pedal to raise 5,000 cc per hook (total 40 lbs) to a height of 9’
- Depress the “DOWN” foot pedal to smoothly lower solution containers, making changeovers a snap
- Heavy cast aluminum base provides stability at maximum and minimum capacity
- 2-column lift system holds a total of 20,000 cc irrigation containers (two 5,000 cc bags per column)
- Battery operation eliminates hazardous electric cords on OR floor

<u>Product</u>	<u>Description</u>	<u>Weight</u>	<u>Height</u>
#741315	... Power Lifter® Irrigation Stand	... 57 lbs.	... Max 9’
#741316	... Low Voltage Light Kit



Keeps O.R. Safe & Operable

The battery charger plugs into the access port on the base of the POWER LIFTER® & POWER LIFTER® II.

The POWER LIFTER® & POWER LIFTER® II will not work while the unit is plugged in. We suggest that a recharge schedule be established to prolong the battery life.



Optional low voltage LED light located in between the foot pedals indicates when unit needs recharging.



Preview Power Lifter® II video at:

www.omnimedbeam.com

