

Slides for Professional Model




Test No. 1, SA-APS-1 Allen Picture Test

Line 1 – Binocularly – Tests performance of both eyes together. Lines 2, 3, 4 – Acuity – Right and Left Eye
Calibrated at 20/30, 20/40 and 20/50.



Test No. 5, SA-BRL-4 – Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, 20/100, and 20/200.
Tests can be administered binocularly or with either eye occluded.

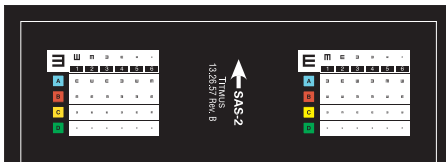


Test No. 2, SAS-2 – Acuity Right and Left Eye

Conventional Tumbling “E” symbol.
Child states which way the “E” points in each square.
Acuity levels of 20/20, 20/30, 20/40 and 20/50.



Test No. 6, SLP-1 – Lateral Phoria

Measures Esophoria or Exophoria in
1 prism diopter steps.

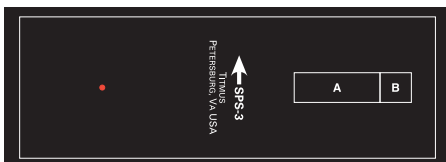

Test No. 3, SPS-3 – Heterophorias

Muscle imbalance – vertical and lateral – is tested at predetermined pass/fail levels.
This test is administered first at far, then at near.



Test No. 7, SVP-1 – Vertical Phoria

Measures right or left Hyperphorias in one-half prism diopter steps.



Test No. 4, SCD-1 – Color Perception

Consists of six accurately and authentically reproduced Ishihara Pseudo-Isochromatic plates.
Detects color deficiency but does not classify.



Test No. 8, SA-BRL-2A – Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, and 20/100.
Tests can be administered binocularly or with either eye occluded.