


Slides for School Model



Test No. 1, SF-1 – Binocularity Both Eyes


To determine whether both eyes are being used together.

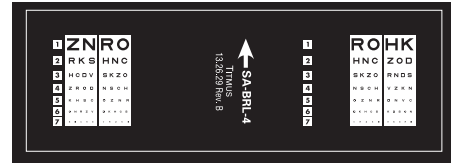

Test No. 5, SA-BRL-4 – Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, 20/100, and 20/200.

Tests can be administered binocularly or with either eye occluded.

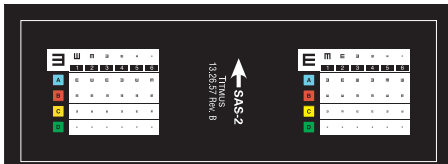


Test No. 2, SAS-2 – Acuity Right and Left Eye

Conventional Tumbling “E” symbol.


Child states which way the “E” points in each square.

Acuity levels of 20/20, 20/30, 20/40 and 20/50.

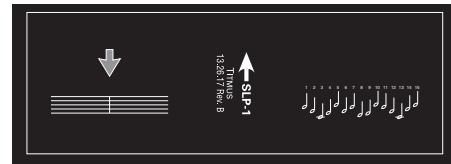

Test No. 6, SLP-1 – Lateral Phoria

Measures Esophoria or Exophoria in 1 prism diopter steps.

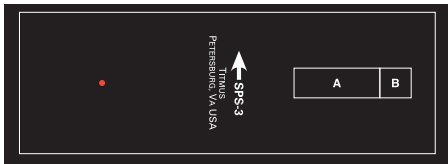

Test No. 3, SPS-3 – Heterophorias


Muscle imbalance – vertical and lateral – is tested at predetermined pass/fail levels.

This test is administered first at far, then at near.



Test No. 7, SVP-1 – Vertical Phoria

Measures right or left Hyperphorias in one-half prism diopter steps.

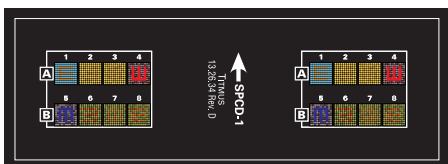

Test No. 4, SPCD-1 – Color Perception

Titmus Pseudo-Isochromatic plates reliably detect gross color anomalies, but do not classify the type of deficiency.



Test No. 8, SA-BRL-2A – Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70 and 20/100.

Tests can be administered binocularly or with either eye occluded.